

***It's Interventional.***

*Patient centred people*

# ***MedOne Ultra automatic biopsy device.***



**BIOPSY**

*hello@itsinterventional.com*

*itsinterventional.com*



# *It's about reliability...*

*MedOne Ultra is a programmable automatic biopsy needle that optimises spring compression for fast and reliable tissue capture even in dense and calcified tissues.*

*Designed with ergonomics and performance in mind, MedOne Ultra combines ease of use with reliable and repeatable sample capture.*

## **Product benefits:**

- **Maximum sample volume capture**
- **Smooth penetration of fibrous and calcified tissues**
- **Lightweight compact design**
- **Single-handed and controlled operation**
- **Precision sampling capability**
- **Suitable for a wide range of procedures**
- **Ergonomically designed handle**
- **Smooth charging**




## **Additional information?**

If you would like further assistance with using MedOne Ultra, or need advice on how to order, please call us on: **0114 268 8880**

**Product features:**

Colour coded loading levers for gauge identification

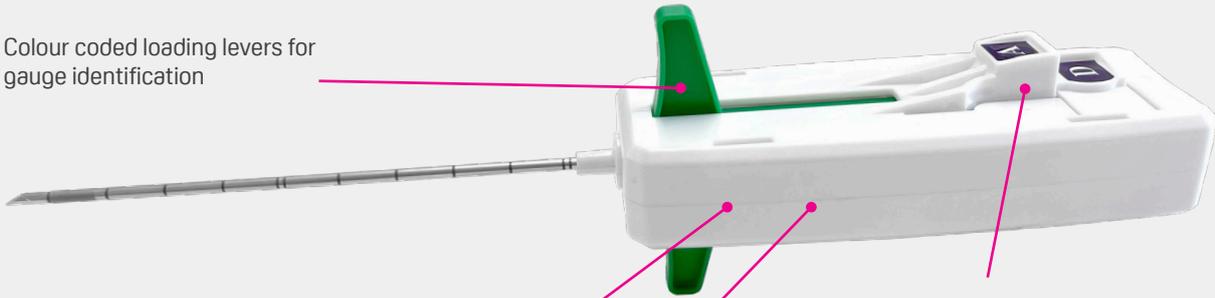


Internal spring firing mechanism provides high spring velocity for tissue penetration and sampling capability

Ergonomic design

Two programmable firing modes

**A** **Automatic** - releases cannula and stylet in rapid sequence and captures tissue sample with the push of a button


**D** **Delay mode** - releases stylet and allows time for tissue to settle in notch; user verifies and documents stylet position, repositions if necessary and activates the cutting cannula

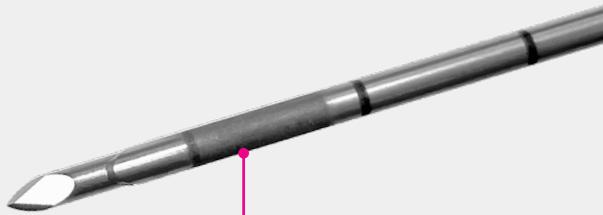

Echogenic treatment for visualisation under ultrasound

Thin-walled sharp cutting cannula

20mm ground sample notch

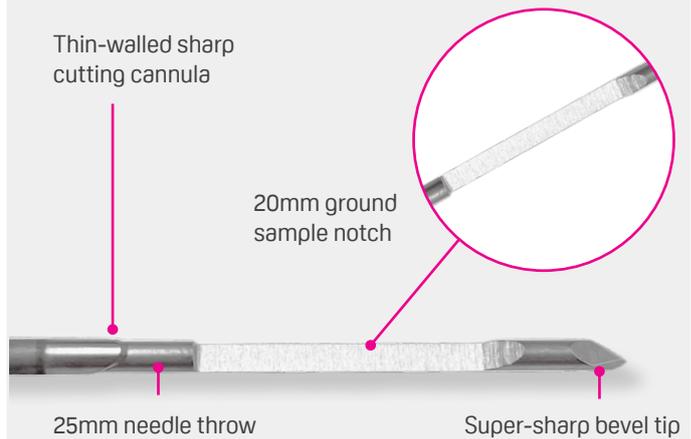

25mm needle throw

Super-sharp bevel tip



Lightweight compact design


Smooth charging action

Single-handed loading for ease of use

**Coaxial Introducer**

Skin stop

cm depth markings

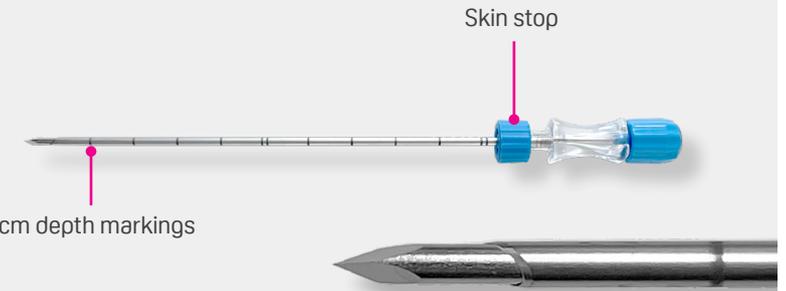

Blunt tip cannula option



## Ordering Information

Product Code	Description	Gauge	Colour	Box Quantity
MU-12-060-00	MedOne Ultra 12g x 6cm	12		5
MU-12-100-00	MedOne Ultra 12g x 10cm			5
MU-12-120-00	MedOne Ultra 12g x 12cm			5
MU-12-160-00	MedOne Ultra 12g x 16cm			5
MU-12-200-00	MedOne Ultra 12g x 20cm			5
MU-14-060-00	MedOne Ultra 14g x 6cm	14		5
MU-14-100-00	MedOne Ultra 14g x 10cm			5
MU-14-120-00	MedOne Ultra 14g x 12cm			5
MU-14-160-00	MedOne Ultra 14g x 16cm			5
MU-14-200-00	MedOne Ultra 14g x 20cm			5
MU-16-060-00	MedOne Ultra 16g x 6cm	16		5
MU-16-100-00	MedOne Ultra 16g x 10cm			5
MU-16-120-00	MedOne Ultra 16g x 12cm			5
MU-16-160-00	MedOne Ultra 16g x 16cm			5
MU-16-200-00	MedOne Ultra 16g x 20cm			5
MU-18-060-00	MedOne Ultra 18g x 6cm	18		5
MU-18-100-00	MedOne Ultra 18g x 10cm			5
MU-18-120-00	MedOne Ultra 18g x 12cm			5
MU-18-160-00	MedOne Ultra 18g x 16cm			5
MU-18-200-00	MedOne Ultra 18g x 20cm			5
MU-18-250-00	MedOne Ultra 18g x 25cm			5
MU-20-060-00	MedOne Ultra 20g x 6cm	20		5
MU-20-100-00	MedOne Ultra 20g x 10cm			5
MU-20-120-00	MedOne Ultra 20g x 12cm			5
MU-20-160-00	MedOne Ultra 20g x 16cm			5
MU-20-200-00	MedOne Ultra 20g x 20cm			5

MedOne Ultra is also available with a packaged coaxial introducer and a packaged coaxial introducer with a blunt tip stylet.


Manufactured by



Distributed by:

It's Interventional  
Albreda House  
Lydgate Lane  
Sheffield, S10 5FH

T: +44 (0) 114 268 8880  
F: +44 (0) 114 268 8881  
orders@itsinterventional.com  
www.itsinterventional.com

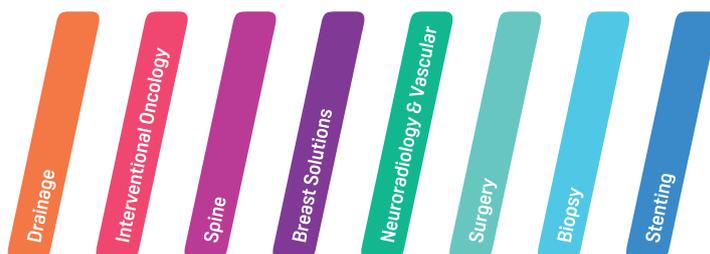

**It's Interventional.**

Patient centred people